

# Targeted Selected Reaction Monitoring Mass Spectrometric Immunoassay for Insulin-like Growth Factor 1

Eric E. Niederkofler<sup>1</sup>, David A. Phillips<sup>1</sup>, Bryan Krastins<sup>2</sup>, Vathany Kulasingam<sup>3</sup>, Urban A. Kiernan<sup>1</sup>, Kemmons A. Tubbs<sup>1</sup>, Scott M. Peterman<sup>2</sup>, Amol Prakash<sup>2</sup>, Eleftherios P. Diamandis<sup>3</sup>, Mary F. Lopez<sup>2</sup>, Dobrin Nedelkov<sup>1\*</sup>

**1** Thermo Fisher Scientific, Tempe, Arizona, United States of America, **2** Thermo Fisher Scientific, BRIMS, Cambridge, Massachusetts, United States of America, **3** University Health Network and University of Toronto, Toronto, Ontario, United States of America

## Abstract

Insulin-like growth factor 1 (IGF1) is an important biomarker of human growth disorders that is routinely analyzed in clinical laboratories. Mass spectrometry-based workflows offer a viable alternative to standard IGF1 immunoassays, which utilize various pre-analytical preparation strategies. In this work we developed an assay that incorporates a novel sample preparation method for dissociating IGF1 from its binding proteins. The workflow also includes an immunoaffinity step using antibody-derivatized pipette tips, followed by elution, trypsin digestion, and LC-MS/MS separation and detection of the signature peptides in a selected reaction monitoring (SRM) mode. The resulting quantitative mass spectrometric immunoassay (MSIA) exhibited good linearity in the range of 1 to 1,500 ng/mL IGF1, intra- and inter-assay precision with CVs of less than 10%, and lowest limits of detection of 1 ng/mL. The linearity and recovery characteristics of the assay were also established, and the new method compared to a commercially available immunoassay using a large cohort of human serum samples. The IGF1 SRM MSIA is well suited for use in clinical laboratories.

**Citation:** Niederkofler EE, Phillips DA, Krastins B, Kulasingam V, Kiernan UA, et al. (2013) Targeted Selected Reaction Monitoring Mass Spectrometric Immunoassay for Insulin-like Growth Factor 1. *PLoS ONE* 8(11): e81125. doi:10.1371/journal.pone.0081125

**Editor:** John Matthew Koomen, Moffitt Cancer Center, United States of America

**Received:** March 26, 2013; **Accepted:** October 9, 2013; **Published:** November 21, 2013

**Copyright:** © 2013 Niederkofler et al. This is an open-access article distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited.

**Funding:** The authors have no support or funding to report.

**Competing Interests:** EEN, DAP, BK, UAK, KAT, SMP, AP, MFL, and DN are employed by Thermo Fisher Scientific. This does not alter the authors' adherence to all the PLOS ONE policies on sharing data and materials.

\* E-mail: dnedelkov@populationproteomics.org

## Introduction

Insulin-like growth factor 1 (IGF1) is an important biological protein implicated in many physiological processes - from cell proliferation, differentiation and apoptosis, to tissue growth and organ specific functions [1]. IGF1 has been traditionally assayed from human plasma or serum samples using enzymatic or radioimmunoassays, but not without some issues - primarily associated with the dissociation from the IGF-binding proteins (IGFBPs) [2], and the dynamic range of the assays [3,4]. The proliferation of mass spectrometry (MS) in the late 90's prompted several efforts toward its implementation in IGF1 detection. In 2001, two groups reported on the development of HPLC electrospray mass spectrometry methods for quantification of IGF1 from standard samples [5,6]. In 2004 Nelson *et al.*, described an immunoaffinity-based MALDI-TOF MS method for quantification of IGF1 from human plasma samples [7]. A few years later, Kirsch *et al.*, developed an IGF1 quantification method using liquid chromatography coupled to isotope dilution mass spectrometry [8]. Several more MS-based methods were published in 2008 and 2009, some of which included IGF1 immunoaffinity retrieval [9], or acetonitrile plasma depletion step in front of the LC-MS/MS [10]. The last few years have seen the convergence of the IGF1 analytical methods toward solid-phase extraction sample preparation and LC-MS [11–15].

The solid phase extraction method was developed for the first IGF1 immunoassays in the 80s and 90s [16–18]. The method, which involves acidic disruption of the IGF1-IGFBP complex and ethanol precipitation of the IGFBPs - leaving free IGF1 in solution, has changed little since then, and most MS-based workflows have incorporated it without any significant modifications. However, it seems that the high-end capabilities of MS could be better matched to a simpler, high-throughput IGF1 sample preparation method. In this work we have adopted and modified the sample preparation approach described by Nelson *et al.* [7], and combined it with a mass spectrometric immunoassay (MSIA) [19,20] method using automated LC-MS/MS for selected reaction monitoring (SRM) and long R3 IGF1 (LR3-IGF1) as an internal reference standard. The resulting IGF1 SRM MSIA was characterized, validated, and employed in screening of clinical human plasma samples, and the data compared to that obtained with a clinical immuno analyzer.

## Materials and Methods

### Reagents

Polyclonal rabbit anti-human IGF1 affinity purified antibody (Cat. No. PA0362), recombinant human IGF1 (Cat. No. C-RI500c), and recombinant human LR3-IGF1 (Cat. No. LRM001) were obtained from Cell Sciences (Canton, MA). Custom MSIA

Pipette Tips (Cat. No. 991CUS02), and phosphate buffered saline (PBS, 28374) were obtained from Thermo Fisher Scientific (San Diego, CA). Ammonium bicarbonate (Cat. No. 09830), calcium chloride (C1016), dithiothreitol (43815), glucagon (G2044), iodoacetamide (I1149), and TWEEN 20 (P7949) were obtained from Sigma-Aldrich (St. Louis, MO). Bovine serum albumin (BSA, Cat. No. 126609), and formic acid (11670) were obtained from EMD Millipore (Billerica, MA). Sterile water (Cat. No. AB02120), acetonitrile (AB00120), isopropyl alcohol (AB00866), SDS (AB01920), and trifluoroacetic acid (TFA, AB02010) were purchased from American Bioanalytical (Natick, MA). Sequencing grade modified trypsin (Cat. No. V511) was obtained from Promega (Madison, WI).

### Human serum and plasma samples

For the method development and validation, several human plasma and serum samples were purchased from ProMedDX (Norton, MA, USA). The samples were received labeled only with a barcode and supplied with an accompanying specification sheet containing information about the gender and age of the donor. For the method comparison studies, 289 serum samples were obtained from the University Health Network (Toronto, ON, Canada), under the approval of the University Health Network Research Ethics Board (REB # 09-0486-T) for the research use of left over routine clinical specimens in the Core Diagnostic Laboratory. Written informed consent from the participants was not required based on the Tri-Council Policy Statement: Ethical Conduct for Research Involving Humans (TCPS2, Article 3.7). The samples were received numbered and without any identifiers, and with their IGF1 levels determined using the Siemens Immulite 2000 IGF-1 assay. All samples were aliquoted and stored at  $-80^{\circ}\text{C}$  until used.

### Preparation of standards and analytical samples

For the standard curve generation, the IGF1 stock (1 g/L) was serially diluted to 1,500; 1,000; 500; 100; 25; 10; 5; and 1 ng/mL, with 10 mM PBS containing 1% (w/v) BSA (standards buffer). The internal reference standard (LR3-IGF1, 1 g/L) was also serially diluted in standards buffer to a final concentration of 500 ng/mL. The analytical samples were prepared in individual wells of a microtiter plate by combining 100  $\mu\text{L}$  sample buffer (10 mM PBS w/0.3% (w/v) SDS), 20  $\mu\text{L}$  of the 500 ng/mL LR3-IGF1 solution, and 40  $\mu\text{L}$  of either the IGF1 standards (for the standard curve generation) or undiluted human plasma or serum. The microplates were then shaken at room temperature for 1 h on an orbital shaker, at 1,000 rpm, to dissociate the IGF1 from the IGF1-IGFBPs.

### Mass spectrometric immunoassay

The immunoaffinity retrieval of IGF1 and LR3-IGF1 from samples was performed using MSIA-Tips derivatized with the IGF1 antibody, which were prepared as previously described [21]. The MSIA-Tips were mounted onto the head of the Versette Automated Liquid Handler (Thermo Fisher Scientific, Hudson, NH) and initially rinsed with assay buffer (10 mM PBS w/0.1% TWEEN 20), with 10 cycles (1 cycle consisting of a single aspiration and dispense of a 100  $\mu\text{L}$  volume,  $\sim 3$  s), from a single 150  $\mu\text{L}$  buffer aliquot placed in the well of a microplate. Next, the MSIA-Tips were immersed into the wells of the microplate containing the samples, and 100 aspirations and dispense cycles were performed (100  $\mu\text{L}$  volumes each), allowing for simultaneous affinity capture of IGF1 and LR3-IGF1. The MSIA-Tips were then rinsed with assay buffer (100 cycles) from another microplate, and twice with water (10 cycles each) from two more microplates (100  $\mu\text{L}$  volumes aspiration and dispenses, from 150  $\mu\text{L}$  placed in

each well). The captured proteins were eluted with 33% acetonitrile/0.4% (v/v) TFA by aspirating and dispensing a 20  $\mu\text{L}$  volume 30 times, from a total of 50  $\mu\text{L}$  in the wells of a microplate. The proteins-containing eluates were dried down in a SpeedVac concentrator until dry, and re-suspended in 30  $\mu\text{L}$  30% N-propanol/100 mM ammonium bicarbonate pH 8.5/10 mM DTT. The microplate was sealed and vortexed for 30 s, spun down, and incubated at  $37^{\circ}\text{C}$  for 30 min to reduce disulfides. A 2.4  $\mu\text{L}$  aliquot of 0.5 M iodoacetamide (in 100 mM ammonium bicarbonate, pH 8.5) was added into each well, and the microplate was sealed again, vortexed for 30 s, spun down, and incubated in the dark at room temperature for 30 min. A 92.5  $\mu\text{L}$  aliquot of 5 mM  $\text{CaCl}_2$ /0.1 M ammonium bicarbonate pH 8 warmed to  $50^{\circ}\text{C}$  was added to each well, and digestion of the reduced and alkylated proteins was initiated by adding 25  $\mu\text{L}$  of 4 mg/L warm ( $50^{\circ}\text{C}$ ) trypsin to each well. The plate was sealed, vortexed for 60 s, spun down and incubated at  $50^{\circ}\text{C}$  for 2 h. The plate was then cooled down at  $4^{\circ}\text{C}$  for 5 min, and 2.3  $\mu\text{L}$  of glucagon (1 g/L in 0.2% formic acid; as a peptide carrier) and 3  $\mu\text{L}$  of 100% formic acid were added to each well. The microplate was positioned in an auto sampler (PAL HTC Accela1-TMO, Thermo Scientific) and the entirety of each sample (155  $\mu\text{L}$ ) was injected onto a 2.1-mm $\times$ 100-mm Accucore aQ 2.6- $\mu\text{m}$  particle C18 column (Thermo Fisher Scientific). Reversed-phase separations were carried out on an Accela 1250 pump (Thermo Fisher Scientific) at a flow rate of 240  $\mu\text{L}/\text{min}$ . Solvent A was 0.2% formic acid in LC-MS grade water, and solvent B was 0.2% formic acid in LC-MS grade acetonitrile (Optima grade reagents, Thermo Fisher Scientific). SRM assays were developed on a TSQ Vantage triple quadrupole mass spectrometer equipped with a HESI-II source (Thermo Fisher Scientific). A mass window of full-width at half maximum of 0.7 (unit resolution) was used in the SRM assays, as a result of the immuno-enriched samples having high signal-to-noise. Pinpoint software (Thermo Fisher Scientific) was used for targeted protein quantification. SRM transitions from the simultaneously immuno-enriched and digested IGF1 (residues 1–21) and LR3-IGF1 (residues 17–34) were monitored, with the peak area ratio of IGF1 to LR3-IGF1 used for quantification.

## Results and Discussion

### Sample preparation

The most common IGF1 sample preparation protocol utilizes acidic conditions to dissociate the IGF1-IGFBPs complex, and ethanol to precipitate IGFBPs from the sample, leaving IGF1 free in solution for subsequent analysis [16–18]. An alternative approach is to release the IGF1 from its complex, and assay the protein directly, without further sample manipulation which might result in some IGF1 loss. Nelson *et al.*, [7] devised an IGF1 sample preparation method that utilized sodium dodecyl sulphate (SDS) to dissociate the IGF1 from the complex with the IGFBPs. We adopted the SDS-based method in the current work, but with a few modifications. First, the optimal SDS concentration at which all of the IGF1 is dissociated from the IGFBPs was re-examined and determined empirically. Next, the SDS concentration was optimized so that it does not interfere with the antibody-antigen interaction (i.e., does not denature the protein or the antibody). This was particularly important because, to avoid any re-association of IGF1 with the IGFBPs, we eliminated any addition of extra assay buffer to the analytical samples (which was part of the method described by Nelson *et al.*, [7]). The final protocol consisted of adding IGF1 sample and LR3-IGF1 standard aliquots to a microplate well containing PBS sample buffer with 0.3% SDS, and shaking the plate for 60 min at room temperature, after which

the IGF1 and LR3-IGF1 were immunoaffinity retrieved from the samples. In this sample preparation method the transfer of reagents is minimized, and the IGF1 sample never leaves the microplate, thus ensuring more accurate IGF1 measurement.

### Internal Reference Standard

Another key feature of the sample preparation is addition of the internal reference standard (IRS) at the beginning of the assay. It is important that the IRS goes through the same processing as does the protein target that is being assayed, to control for possible losses during these processes [22]. In effect, when properly designed and added at the beginning of the assay, the IRS serves as a normalizer for all sample processing and data acquisition steps – analyte extraction, reduction/alkylation, digestion and MS response. In contrast, most targeted proteomics assays advocate the use of surrogate peptides and their corresponding isotope labeled counterparts as internal reference standards [23,24], which are added to the sample *after* the initial tryptic digest of the plasma, and can lead to incorrect quantitative results due to incomplete digestion of the targeted protein. The use of isotope-labeled protein standards for quantification (added to the sample prior to processing) has been suggested as a better alternative [25], but they remain expensive and difficult to manufacture. In contrast, readily available protein analogues can also serve as adequate IRS. We have used such analogues in the past, and most of them have been proteins from other species that differ very little in the amino acid sequence [7,26,27]. In this work we utilized long arginine 3-IGF1 (LR3-IGF1) as an internal reference standard. LR3-IGF1 differs from native IGF1 in that it has an arginine instead of glutamic acid at position 3 in the amino acid sequence, and it contains extra 13 amino acids at the N-terminus (MFPAMPLLSLNFVN), the net result of which is reduced binding of LR3-IGF1 to the IGF1BPs. Other possibilities for an IRS include rat IGF1 [13], mouse IGF1, and R3-IGF1 (which will result in identical transition pairs in the MS/MS as those of LR3-IGF1), all of which are readily available as recombinant proteins. We determined empirically that the sequence modifications in LR3-IGF1 did not influence the binding to the IGF1 antibody that was used in the subsequent steps of immunoaffinity purification, which is an important prerequisite.

### Assay Development and Optimization

We opted for an immunoaffinity approach to IGF1 preparation for the subsequent MS detection. Medium- to high-concentration

plasma proteins (>1 mg/L) can be detected via LC-MS/MS without enrichment [28], but for most others an immuno-enrichment is more often than not a necessary step. Inter-laboratory evaluation of immunoaffinity enrichment peptide multiple reaction monitoring (MRM) MS assays has demonstrated limits of detection of 1 ng/mL [23]. In our past work we have shown that with the use of appropriate devices for immunoaffinity purification, lowest limits of detection (LLOD) of 10 pg/mL or less are possible [19,29]. The immunoaffinity capture of IGF1 and LR3-IGF1 was achieved through the use of MSIA-Tips – disposable automation research tips fitted with a small porous microcolumn at the distal end, onto which an IGF1 antibody was immobilized using standard protocols [21]. The optimal amount of antibody was empirically determined to be 1 µg/tip, considering both the performance of the assay (primarily the LLOD) and the cost factor for the antibody. The IGF1 affinity capture and elution protocol were designed with simplicity in mind, containing a minimum number of steps needed for successful capture of the proteins from the sample. When performed on the Versette Automated Liquid Handler, the affinity retrieval and elution of IGF1 and LR3-IGF1 was executed in <45 min, for 96 samples at a time. The subsequent trypsin digestion and LC/MS protocols were adapted from our previous work [19], with few modifications as described in the Methods section. The MS/MS transitions that were monitored for the quantification of IGF1 and LR3-IGF1 are listed in **Table 1**.


### Assay Performance

Standard curves were generated by plotting the IGF1/LR3-IGF1 peak area ratios vs. the concentration of the IGF1 standards. A typical IGF1 standard curve (**Fig. 1**) covered the range from 1 to 1,500 ng/mL IGF1, with good linearity ( $R^2 = 0.999$ ,  $SEE = 0.135$ ). The lowest limit of detection (LLOD) achieved with the assay was 1 ng/mL (equating to ~5 femtomoles of IGF1), which was determined as 2 standard deviations above background noise from blank samples. A lower limit of quantifications (LLOQ) of 5 ng/mL was also achieved, as the lowest concentration that was reproducibly quantified. The intra-assay precision (within-run) was determined by analyzing three plasma samples, in triplicates, each with a single standard curve. The inter-assay precision (run-to-run) was determined by analyzing one plasma sample three times, on different days, with separate standard curves each time. The results are shown in **Table 2**, and indicate CVs of less than

**Table 1.** MS/MS transitions for IGF1 and LR3-IGF1 used for the quantification.

IGF1	GPETLC(carboxymethyl)GAELVDALQFVC(carboxymethyl)GDR	
	770.32 (+3)	508.181
		607.249
		754.318
		882.377
		995.461
LR3-IGF1	TLC(carboxymethyl)GAELVDALQFVC(carboxymethyl)GDR	
	675.98 (+3)	347.167
		508.181
		607.249
		754.318
		882.377

doi:10.1371/journal.pone.0081125.t001

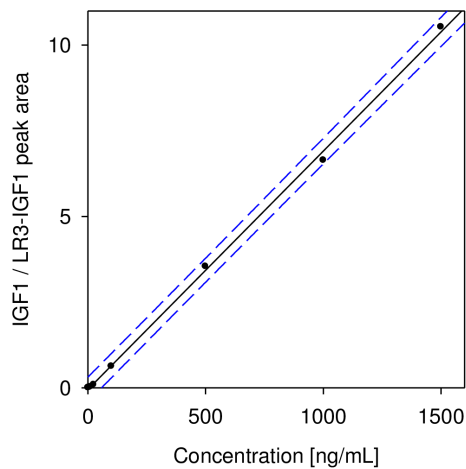

**Figure 1. Representative standard curve for the IGF1 SRM-MSIA.** Plotted are the peak area ratios of IGF1/LR3-IGF1 over the standards concentrations. Solid line: linear fit with  $R^2=0.999$ ,  $SEE=0.135$ . Dotted lines: 95% prediction intervals. doi:10.1371/journal.pone.0081125.g001

10%. To determine the linearity of the assay, plasma samples with known IGF1 concentrations were serially diluted, analyzed with the mass spectrometric immunoassay to determine the IGF1 concentrations, and the results compared to those expected (**Table 3**). Spiking recovery experiments were also performed by spiking plasma samples, with known low IGF1 concentrations, with increasing amounts of recombinant IGF1, followed by analysis with the assay to determine the total IGF1 concentration, and comparison of the results with those expected (**Table 4**).

### Method Comparison

The fully-developed, characterized, and validated IGF1 SRM-MSIA assay was then used to analyze a set ( $n=289$ ) of human serum samples. The samples were provided by the University Health Network, with their IGF1 levels determined by the Siemens Immulite 2000 Platform [30]. The IGF1 concentrations determined with the SRM-MSIA correlated well with those obtained with the Immulite System, with a Passing Bablock fit [31] of  $-9.11+0.82x$ . The Bland-Altman plot [32] shows a slight negative bias of 26% (**Fig. 2**). This is not unusual, and it has been observed in a number of method comparison studies [12,33]. Further studies might be needed to delineate the source of these deviations, which is particularly important for samples with high IGF1 (e.g. in patients with growth hormone excess - acromegaly; similarly, the effects of increased IGFBP levels on the performance of the assay would need to be further evaluated). The IGF1 SRM-MSIA assay LLOD (1 ng/mL) is lower than the LOD of the Immulite assay (20 ng/mL), and that of a recently published LC-

**Table 2. Intra-and inter-assay precision.**

Intra-assay CVs				Inter-assay CV	
Sample	1	2	3	STDEVP	11.1
STDEVP:	6.68	9.93	5.88	MEAN (ng/mL):	151
MEAN (ng/mL):	166	149	139	CV (%)	7.36
CV (%):	4.03	6.67	4.24		

doi:10.1371/journal.pone.0081125.t002

**Table 3. Assay linearity.**

Sample	Dilution	Observed	Expected	Recovery
		ng/mL	ng/mL	O/E%
1		306		
	2×	146	153	95.4
	4×	72.7	77.0	94.4
	8×	39.1	38.0	103
2		239		
	2×	131	119	110.1
	4×	57.0	60.0	95.0
	8×	34.0	30.0	113

doi:10.1371/journal.pone.0081125.t003

MS/MS assay by Bystrom *et al.* -4 ng/mL [12]. Before the IGF1 SRM-MSIA assay can find clinical application it will need calibrating to the latest IGF-1 international standard (02/254) [34].


### Conclusions

Mass spectrometry-based workflows offer a viable alternative to standard IGF1 immunoassays. While some MS-based methods have been developed for assaying IGF1 from blood spots [35] or urine [15], human plasma remains the most clinically important sample for IGF1 analysis, which requires pre-processing steps in order to contend with IGFBPs. In this work, we have taken the approach of disrupting IGF1/IGFBP complexes *in situ*, as well as streamlined post-processing steps in order to create a highly accurate and reproducible IGF1 SRM MSIA. Critical to the assay is the inclusion of the IRS prior to any pre-processing, extraction, post-processing and analysis steps. Thus, with a single, judiciously chosen IRS, all sources of assay variability can be normalized. Importantly, we have utilized immunoaffinity purification in a *standardized* format to enhance the performance of the assay [36]. While one might argue that immuno-enrichment is not necessary with regards to IGF1, we found it to be necessary for the optimized sample preparation workflow presented in this work (the influence of heterophilic antibodies can be eliminated via the use of blocking tubes [37]). Furthermore, immunoaffinity purification is a prerequisite for detection of most low-abundance proteins, as well as variants of those high- and medium-abundance proteins that

**Table 4. Spiking recovery.**

Sample	Observed	Expected	Recovery
	ng/mL	ng/mL	O/E%
1	117		
	229	217	106
	385	317	121
	588	517	114
2	107		
	322	294	110
	546	482	113
	757	857	88.3

doi:10.1371/journal.pone.0081125.t004

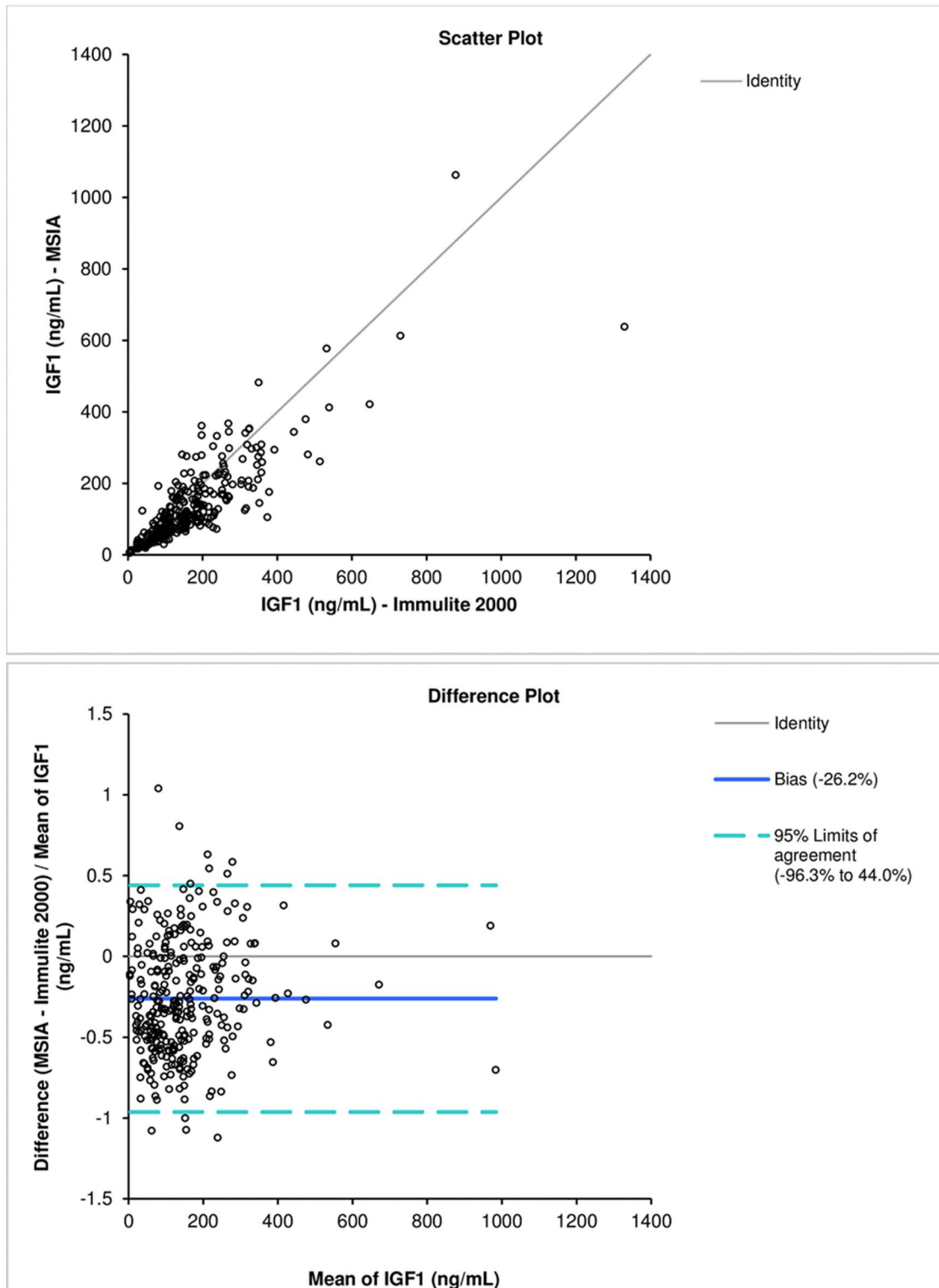

**Figure 2. IGF1 methods comparison.** a) Scatter plot; b) Difference plot.  
doi:10.1371/journal.pone.0081125.g002

are typically found at much lower concentration than their wild-type counterpart. Mass spectrometry, when coupled with immunoaffinity separation, has the necessary sensitivity, automation, and high-throughput capability needed for running hundreds of samples per day. Lastly, in this work we performed trypsin

digestion (post-IGF1 capture and elution) only because the triple quadrupole mass spectrometer employed has a better sensitivity at the lower peptide mass range. Others have shown that digestion is not needed for a MS-based IGF1 assay detection [7,13], and we certainly recognize the fact that without digestion the entire IGF1

assay is much simpler. In that respect we expect to adopt the sample preparation and immunoaffinity steps described in this work to a workflow and assay that employs an advanced MS platform with matched sensitivities to detect intact IGF1. Such assays are certain to find use in the clinical laboratories.

## References

- Zapf J, Froesch ER, Schmid C (1999) Metabolic effects of IGFs. In: Conn M, editor. *Contemporary Endocrinology 17: the IGF system - molecular biology, physiology and clinical applications*. New Jersey: Humana Press.
- Forbes BE, McCarthy P, Norton RS (2012) Insulin-like growth factor binding proteins: a structural perspective. *Front Endocrinol (Lausanne)* 3: 38.
- Clemmons DR (2007) IGF-I assays: current assay methodologies and their limitations. *Pituitary* 10: 121–128.
- Frystyk J, Freda P, Clemmons DR (2010) The current status of IGF-I assays—a 2009 update. *Growth Horm IGF Res* 20: 8–18.
- Bobin S, Popot MA, Bonnaire Y, Tabet JC (2001) Approach to the determination of insulin-like-growth-factor-I (IGF-I) concentration in plasma by high-performance liquid chromatography-ion trap mass spectrometry: use of a deconvolution algorithm for the quantification of multiprotonated molecules in electrospray ionization. *Analyst* 126: 1996–2001.
- de Kock SS, Rodgers JP, Swanepoel BC (2001) Growth hormone abuse in the horse: preliminary assessment of a mass spectrometric procedure for IGF-1 identification and quantitation. *Rapid Commun Mass Spectrom* 15: 1191–1197.
- Nelson RW, Nedelkov D, Tubbs KA, Kiernan UA (2004) Quantitative mass spectrometric immunoassay of insulin like growth factor 1. *J Proteome Res* 3: 851–855.
- Kirsch S, Widart J, Louette J, Focant JF, De Pauw E (2007) Development of an absolute quantification method targeting growth hormone biomarkers using liquid chromatography coupled to isotope dilution mass spectrometry. *J Chromatogr A* 1153: 300–306.
- Bredheft M, Schanzer W, Thevis M (2008) Quantification of human insulin-like growth factor-1 and qualitative detection of its analogues in plasma using liquid chromatography/electrospray ionisation tandem mass spectrometry. *Rapid Commun Mass Spectrom* 22: 477–485.
- Kay RG, Barton C, Velloso CP, Brown PR, Bartlett C, et al. (2009) High-throughput ultra-high-performance liquid chromatography/tandem mass spectrometry quantitation of insulin-like growth factor-I and leucine-rich alpha-2-glycoprotein in serum as biomarkers of recombinant human growth hormone administration. *Rapid Commun Mass Spectrom* 23: 3173–3182.
- Barton C, Kay RG, Gentzer W, Vitthum F, Pleasance S (2010) Development of high-throughput chemical extraction techniques and quantitative HPLC-MS/MS (SRM) assays for clinically relevant plasma proteins. *J Proteome Res* 9: 333–340.
- Bystrom C, Sheng S, Zhang K, Caulfield M, Clarke NJ, et al. (2012) Clinical utility of insulin-like growth factor 1 and 2; determination by high resolution mass spectrometry. *PLoS One* 7: e43457.
- Bystrom CE, Sheng S, Clarke NJ (2011) Narrow mass extraction of time-of-flight data for quantitative analysis of proteins: determination of insulin-like growth factor-1. *Anal Chem* 83: 9005–9010.
- Kay R, Halsall DJ, Annamalai AK, Kandasamy N, Taylor K, et al. (2013) A novel mass spectrometry-based method for determining insulin-like growth factor 1: assessment in a cohort of subjects with newly diagnosed acromegaly. *Clin Endocrinol (Oxf)* 78: 424–430.
- Thomas A, Kohler M, Schanzer W, Delahaut P, Thevis M (2011) Determination of IGF-1 and IGF-2, their degradation products and synthetic analogues in urine by LC-MS/MS. *Analyst* 136: 1003–1012.
- Bang P, Eriksson U, Sara V, Wivall IL, Hall K (1991) Comparison of acid ethanol extraction and acid gel filtration prior to IGF-I and IGF-II radioimmunoassays: improvement of determinations in acid ethanol extracts by the use of truncated IGF-I as radioligand. *Acta Endocrinol (Copenh)* 124: 620–629.
- Breier BH, Gallaher BW, Gluckman PD (1991) Radioimmunoassay for insulin-like growth factor-I: solutions to some potential problems and pitfalls. *J Endocrinol* 128: 347–357.
- Daughaday WH, Mariz IK, Blethen SL (1980) Inhibition of access of bound somatomedin to membrane receptor and immunobinding sites: a comparison of radioreceptor and radioimmunoassay of somatomedin in native and acid-ethanol-extracted serum. *J Clin Endocrinol Metab* 51: 781–788.
- Lopez MF, Rezai T, Sarracino DA, Prakash A, Krastins B, et al. (2010) Selected reaction monitoring-mass spectrometric immunoassay responsive to parathyroid hormone and related variants. *Clin Chem* 56: 281–290.
- Nedelkov D (2012) Mass spectrometry-based protein assays for in vitro diagnostic testing. *Expert Rev Mol Diagn* 12: 235–239.
- Trenchevska O, Kamcheva E, Nedelkov D (2010) Mass spectrometric immunoassay for quantitative determination of protein biomarker isoforms. *J Proteome Res* 9: 5969–5973.
- Lange V, Picotti P, Domon B, Aebersold R (2008) Selected reaction monitoring for quantitative proteomics: a tutorial. *Mol Syst Biol* 4: 222.
- Kuhn E, Whiteaker JR, Mani DR, Jackson AM, Zhao L, et al. (2012) Interlaboratory evaluation of automated, multiplexed peptide immunoaffinity enrichment coupled to multiple reaction monitoring mass spectrometry for quantifying proteins in plasma. *Mol Cell Proteomics* 11: M111 013854.
- Whiteaker JR, Zhao L, Abbatiello SE, Burgess M, Kuhn E, et al. (2011) Evaluation of large scale quantitative proteomic assay development using peptide affinity-based mass spectrometry. *Mol Cell Proteomics* 10: M110 005645.
- Brun V, Dupuis A, Adrait A, Marcellin M, Thomas D, et al. (2007) Isotope-labeled protein standards: toward absolute quantitative proteomics. *Mol Cell Proteomics* 6: 2139–2149.
- Niederkofer EE, Tubbs KA, Gruber K, Nedelkov D, Kiernan UA, et al. (2001) Determination of beta-2 microglobulin levels in plasma using a high-throughput mass spectrometric immunoassay system. *Anal Chem* 73: 3294–3299.
- Tubbs KA, Nedelkov D, Nelson RW (2001) Detection and quantification of beta-2-microglobulin using mass spectrometric immunoassay. *Anal Biochem* 289: 26–35.
- Addona TA, Abbatiello SE, Schilling B, Skates SJ, Mani DR, et al. (2009) Multi-site assessment of the precision and reproducibility of multiple reaction monitoring-based measurements of proteins in plasma. *Nat Biotechnol* 27: 633–641.
- Niederkofer EE, Kiernan UA, O'Rear J, Menon S, Saghir S, et al. (2008) Detection of Endogenous B-Type Natriuretic Peptide at Very Low Concentrations in Patients With Heart Failure. *Circulation: Heart Failure* November 1: 258–264.
- Soldin OP, Dahlin JR, Gresham EG, King J, Soldin SJ (2008) IMMULITE 2000 age and sex-specific reference intervals for alpha fetoprotein, homocysteine, insulin, insulin-like growth factor-1, insulin-like growth factor binding protein-3, C-peptide, immunoglobulin E and intact parathyroid hormone. *Clin Biochem* 41: 937–942.
- Passing H, Bablok (1983) A new biometrical procedure for testing the equality of measurements from two different analytical methods. Application of linear regression procedures for method comparison studies in clinical chemistry, Part I. *J Clin Chem Clin Biochem* 21: 709–720.
- Bland JM, Altman DG (1999) Measuring agreement in method comparison studies. *Stat Methods Med Res* 8: 135–160.
- Krebs A, Wallaschofski H, Spilke-Liss E, Kohlmann T, Brabant G, et al. (2008) Five commercially available insulin-like growth factor I (IGF-I) assays in comparison to the former Nichols Advantage IGF-I in a growth hormone treated population. *Clin Chem Lab Med* 46: 1776–1783.
- Burns C, Rigby P, Moore M, Rafferty B (2009) The First International Standard For Insulin-like Growth Factor-1 (IGF-1) for immunoassay: preparation and calibration in an international collaborative study. *Growth Horm IGF Res* 19: 457–462.
- Cox HD, Rampton J, Eichner D (2013) Quantification of insulin-like growth factor-1 in dried blood spots for detection of growth hormone abuse in sport. *Anal Bioanal Chem* 405: 1949–1958.
- Prakash A, Rezai T, Krastins B, Sarracino D, Athanas M, et al. (2012) Interlaboratory reproducibility of selective reaction monitoring assays using multiple upfront analyte enrichment strategies. *J Proteome Res* 11: 3986–3995.
- Brugts MP, Luermans JG, Lentjes EG, van Trooyen-van Vrouwerff NJ, van der Horst FA, et al. (2009) Heterophilic antibodies may be a cause of falsely low total IGF1 levels. *Eur J Endocrinol* 161: 561–565.

## Author Contributions

Conceived and designed the experiments: DN EEN BK. Performed the experiments: EEN DAP BK. Analyzed the data: EEN DAP BK MFL DN. Contributed reagents/materials/analysis tools: VK UAK KAT SMP AP EPD. Wrote the paper: DN.